




Operating Instructions

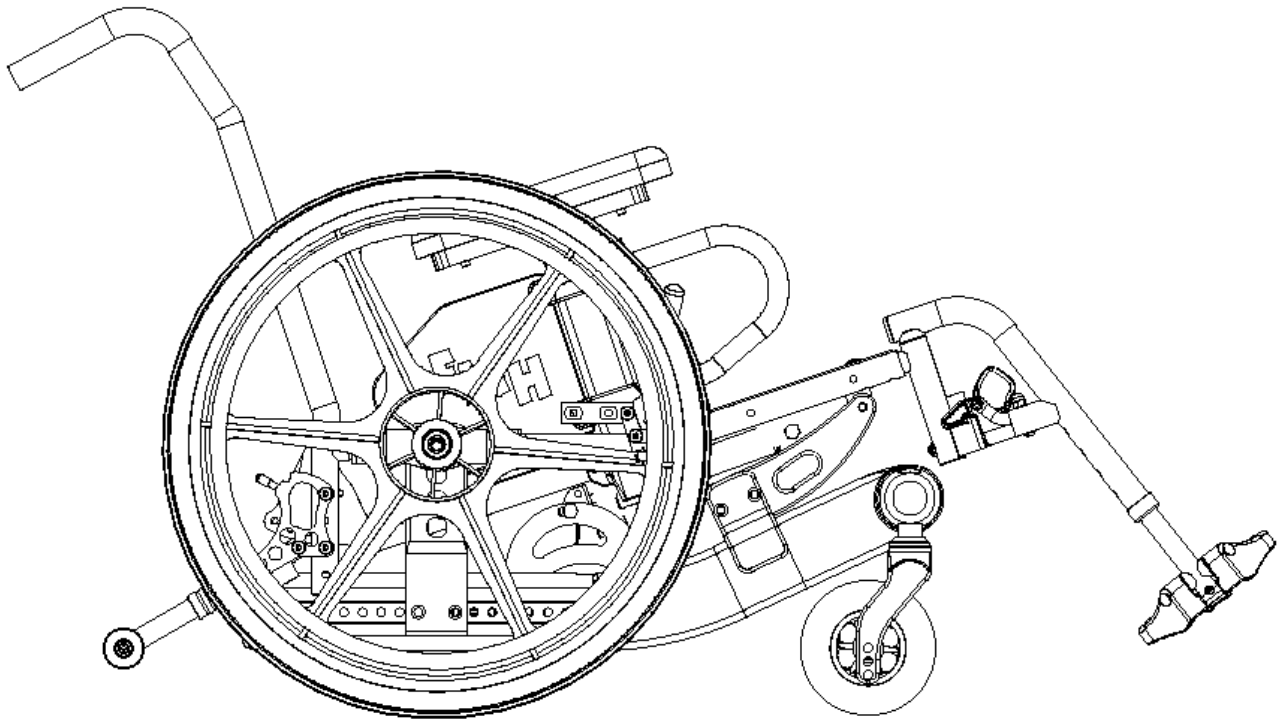

NOTE: Check all parts for shipping damages before using.
In case of damage, do NOT use the equipment.
Contact the Equipment Supplier for further instructions.

WARNING

DO NOT install this equipment without first reading and understanding this instruction booklet. If you are unable to understand these instructions, contact a healthcare professional, dealer or technical personnel before attempting to install this equipment - otherwise, injury or damage may occur.

NOTICE

Information contained within this document is subject to change without notice.

OPERATING SAFETY WARNINGS

Please read and obey all instructions and warnings listed in this manual, ignoring these warnings could result in serious injury to the patient or attendants.

- Ensure anti-tippers are always installed in the correct position when using the wheelchair.
- **Never** lift the wheelchair by the armrests or front rigging.
- Any change of seat depth may require repositioning of the rear wheels, always test the stability of the chair before use.
- **Never** stand on the footplates during transfers in and out of the wheelchair.
- Ensure wheel locks are fully engaged while stopped on any incline/decline and while patient is transferring to and from the wheelchair.
- **Never** tilt the wheelchair on two wheels without assistance. When in any tilt or recline position, this chair must be operated by an attendant.
- **Never** use the wheel locks to stop a moving wheelchair.
- It is not recommended to transport users in vehicles of any kind while seated in the wheelchair. Do not use any kind of wheelchair Tie-Down restraints that are not approved by the FDA.

CAPELLA STANDARD FEATURES

- Tilt Function (17 Degrees)
- 16" – 19" wide
- 15" – 19" depth
- Stroller Bar Low Fixed Back
- Fixed Height Full-Length Arms
 - Height of 9 ½"
 - Lap Tray placement
- 16" STF
- 20" Rear wheels
- 6" Casters
- Adjustable Axle Plate
- 70 Degree Upgraded Dual Swingaway Front Rigging with Footplate
- Spring Back Push to Lock Wheel Locks
- Choice of 8 spirited colors (Black, Blue, Green, Orange, Purple, Red, Silver, Yellow)

* User Weight limit: 350 lbs

SPECIFICATIONS:

Tilt Range: 3 - 20 degrees

Overall Dimensions

(18" x 16" tilt , 16" STF)

Chair Width: 28.75"

Chair Length: 31"

Chair Height: 32"

Chair Weight: 40 Lbs.

(without seat pan, front riggings, wheels, anti tippers and armrests)

SAFETY INSPECTION CHECKLIST

Initial adjustments should be made to suit your personal body structure needs and preference. Thereafter follow these maintenance procedures:

Item	Initially	Weekly	Monthly	Periodically
GENERAL <ul style="list-style-type: none"> Wheelchair rolls straight (no excessive drag or pull to one side) 	X		X	
FRAME AND CROSSING TUBES <ul style="list-style-type: none"> Inspect for loose or missing hardware Inspect for bent frame or crosstubes 	X X			X X
WHEEL LOCKS <ul style="list-style-type: none"> Do not interfere with tires when rolling Pivot points free of wear and looseness Wheel locks easy to engage Wheel locks prevent chair from moving when engaged 	X X X X	X	X X	X
SEAT AND BACK <ul style="list-style-type: none"> Inspect for rips or sagging Inspect for loose or broken hardware Inspect cane and hand grips for wear/looseness 	X X X		X X	X
TIRES <ul style="list-style-type: none"> Inspect for flat spots, cracks and wear Caution: as with any vehicle, the wheels and tires should be checked periodically for cracks and wear and should be replaced. 	X	X		X
REAR WHEELS <ul style="list-style-type: none"> If equipped, quick-release axles lock properly No excessive side movement or binding when lifted and spun Inspect for cracked, bent or broken spokes 	X X X	X X		X
HANDRIMS <ul style="list-style-type: none"> Inspect for signs of rough edges or peeling 	X			X
FRONT CASTERS/FORKS <ul style="list-style-type: none"> Inspect caster fork assembly for proper tension by spinning caster; caster should come to a gradual stop Adjust bearing system if wheel wobbles or binds to a stop. Ensure wheel bearings are clean and free of moisture. Check stem caster journal and lock nut for tightness Inspect casters for cracks and wear Inspect for cracked, bent or broken spokes 	X X X X X	X X	X X	X X
CLEANING <ul style="list-style-type: none"> Clean upholstery and armrests 	X			X

TROUBLESHOOTING

Chair Veers Right/Left	Capeela Wheels	Sluggish Turn or Performance	Caster Flutter	Squeaks and Rattles	Looseness in Chair	Solutions
		X	X	X	X	Check for loose nuts and bolts
X	X		X			Check angle adjustable caster assembly
X	X					Check that rear wheels are equally spaced away from seat frame.

CAPELLA LABELS

WEIGHT CAPACITY
350 LBS. (159 kgs)
 REFER TO OWNER'S MANUAL

WARNING
DO NOT OPERATE WITHOUT
THE ANTI-TIP MECHANISM
IN PLACE

TILT USE CAUTION
WHEN PROPELLING THIS CHAIR
IN THE TILT POSITION. USER
WEIGHT 350 LBS MAXIMUM

OPERATING INSTRUCTIONS

Activating the Tilt Mechanism (refer to Figure 1)

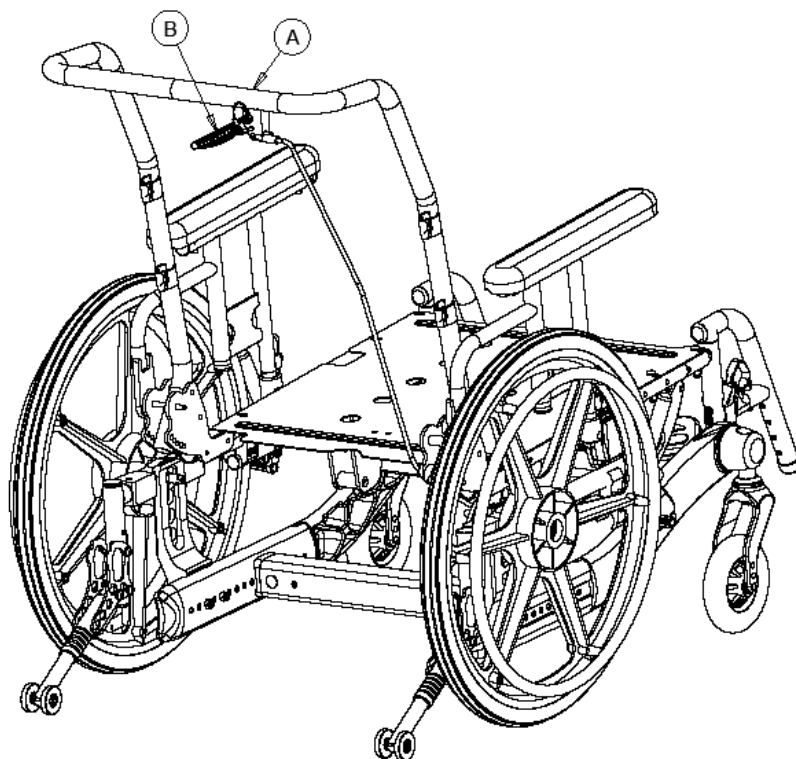


Figure 1

1. Grasp the tilt handle ('B') located on the stroller handle ('A')
2. Pull up on handle to release the wheelchair from the 'locked' position and gently push down on the stroller handle to tilt the chair.
3. When the chair tilts to the desired angle, release the tilt handle and the chair will 'lock' in place at the specific angle.
4. To return the wheelchair to the nominal 'zero' position grasp the tilt handle and pull up on the stroller handle to tilt the chair.
5. Release the tilt handle when the wheelchair reaches the desired angle.

T-Style Fixed Arm Rests (refer to Figure 2a)

Installing and removing T-Style Fixed Armrest

1. Push armrest into the armrest base until lock is engaged into slot ('A')
 2. For removal, push and hold release lever ('B'), while lifting armrest from the armrest base
- Note:** The T-Style Fixed Arm Rest is not adjustable in height

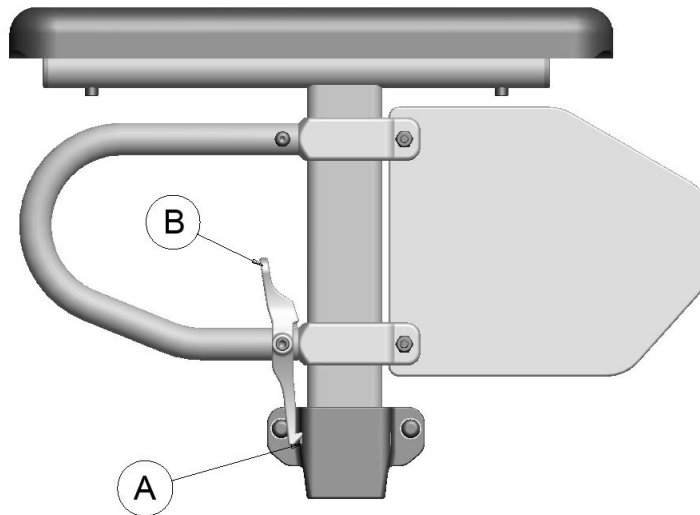


Figure 2a

T-Style Arm Rests – Pin Release(refer to Figure 2b)

ADJUSTING T-STYLE ARM HEIGHT

1. Swing out the height adjustment latch lever to disengage the holding pin.
2. Lift armrest until desired position is acquired (ranges from 9"– 13" in ½" increments).
3. Swing in the height adjustment latch lever to engage the holding pin.
4. Repeat this procedure for the opposite T-style arm assembly.

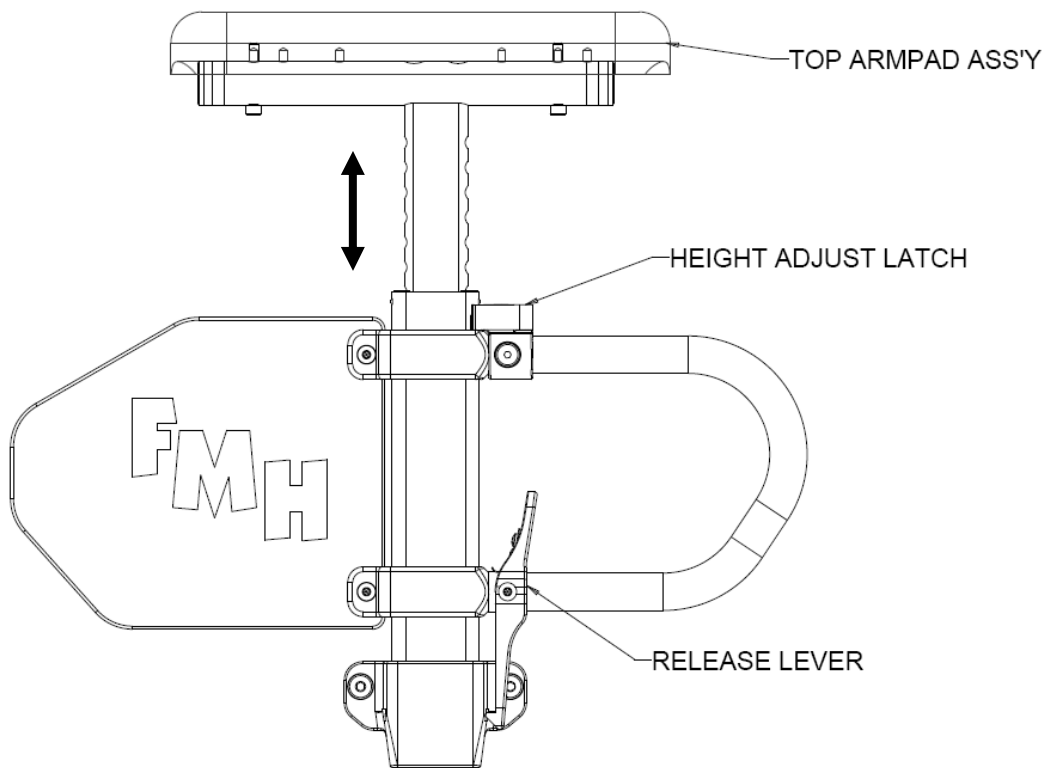


Figure 2b

Flip-up Arm Rests (refer to Figure 2c)

Engaging and Disengaging Flip-up style armrest

1. To engage, push armrest into arm socket until release lever ('A') is engaged into slot.
2. For flip-up, push and hold release lever ('A'), while lifting armrest from arm socket.

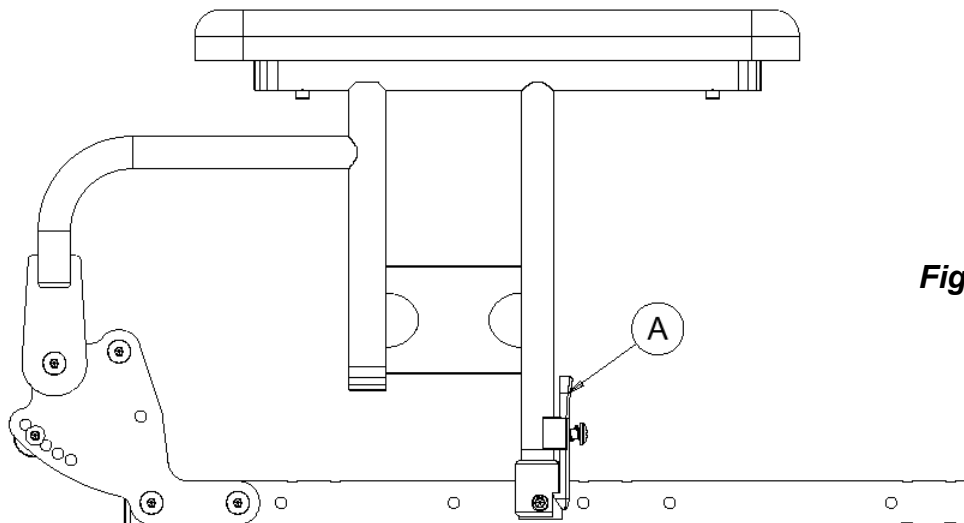


Figure 2c

Footrest Upgraded Dual Swingaway(refer to Figure 3a)

1. Position swing-away footrest onto hanger ('A') ensuring pivot plug sits flush onto hanger tube.
2. To remove footrest pull up handle ('B') releasing latch ('C') and swing footrest outward.

Footplate Length Adjustment

1. Loosen footrest bracket bolt from hole 'D'. Slide footrest extension tube up or down to desired height.
2. Tighten bolt.

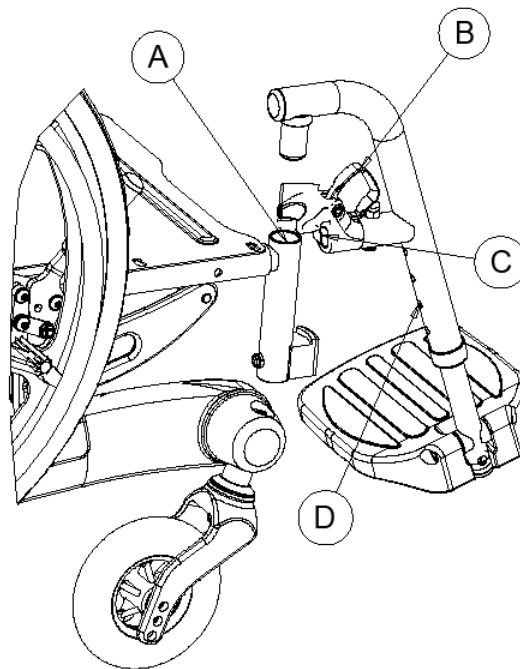


Figure 3a

Elevating Footrest Upgraded Dual Swingaway(refer to Figure 3b)

1. Position swing-away footrest onto hanger ('A') ensuring pivot plug sits flush onto hanger tube.
2. To remove footrest pull up handle ('B') releasing latch ('C') and swing footrest outward.

Footplate Length Adjustment

1. Loosen footrest bracket bolt from slot ('D'). Slide footrest extension tube up or down to desired height.
2. Tighten bolt.

Changing Angle of Elevating Legrest

1. Turn knob ('E') counterclockwise to release the elevating legrest from being locked at an angle.
2. Grasp the elevating legrest outer part and slowly pull towards oneself. The angle of the elevating legrest outer part should become larger.
3. Set at desired angle and then turn knob ('E') clockwise to 'lock' the outer elevating legrest part in place.
4. Repeat for the other side if desired. Ensure that the legrest is securely in place once adjustments have been made.

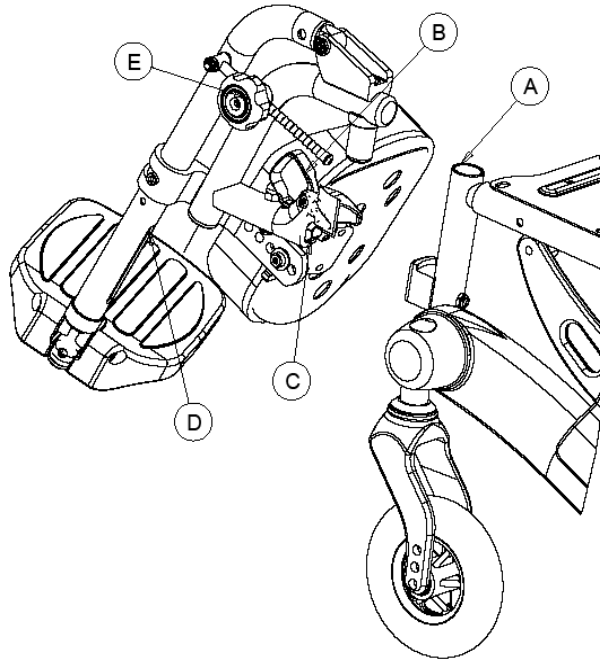


Figure 3b

Chair Back

Fixed Back Option

To adjust Fixed Back angle (See Figure 4a)

Manufacturer's setting is 90 degrees, with each hole representing 7 degree increments (holes range from 83 degrees to 111 degrees).

1. Remove nut and bolt on both sides of the back post.
2. Move back cane to desired angle and replace nuts and bolts making sure covered spacers are in place.

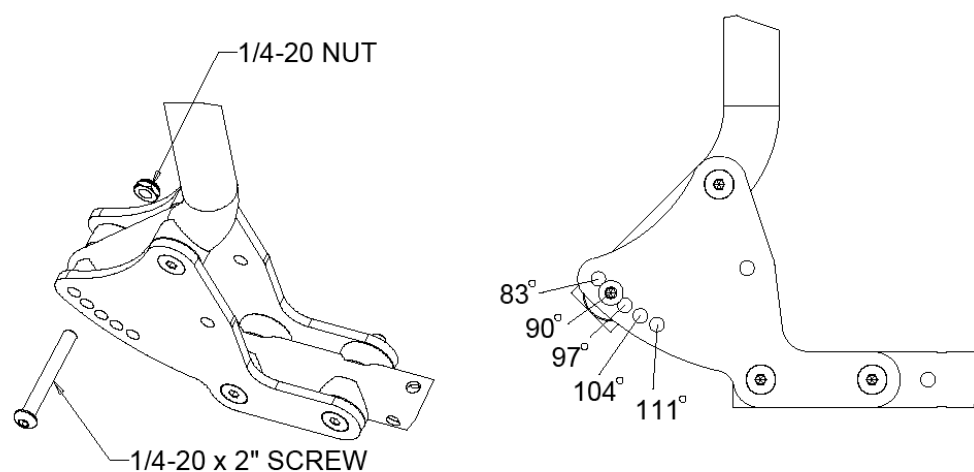


Figure 4a

Fold Down Back Option

To adjust Fold Down Back angle (See Figure 4b)

Manufacturer's setting is 90 degrees, with each hole representing 7 degree increments (holes range from 83 degrees to 111 degrees).

1. Push up on Locking Pin and fold the back away from the holes located for setting the back angle.
2. Remove the shoulder bolt and nut on bolt sides of the back post
3. Place the shoulder bolt and nut in the desired hole location as per the customer specifications. One can refer to the above diagram for hole placement from 83° to 111°.
4. Push up on Locking Pin and fold the back in the new hole locations. Check the angle of the back and adjust the hole location if needed. Ensure the Locking device is engaged and the back is not loose.

To Fold Down the Back

One can fold down the back for easy transport or other reasons by simply pushing up on the locking pin as shown in the diagram below.

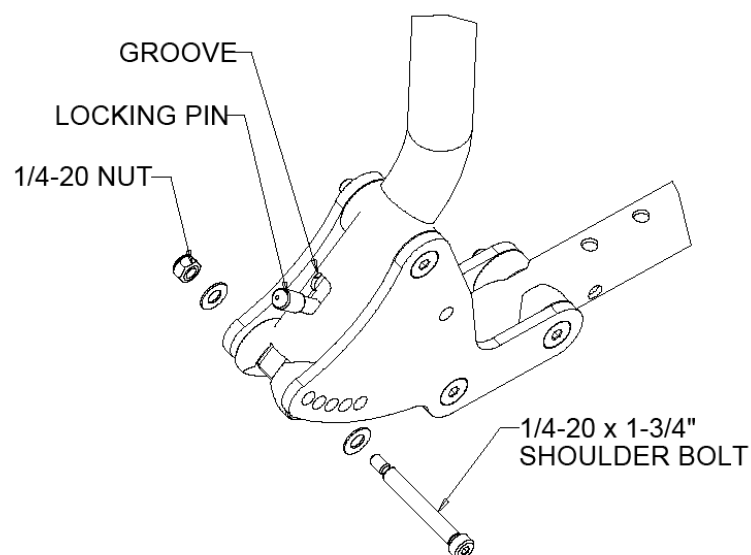


Figure 4b

Changing Seat Width and Depth

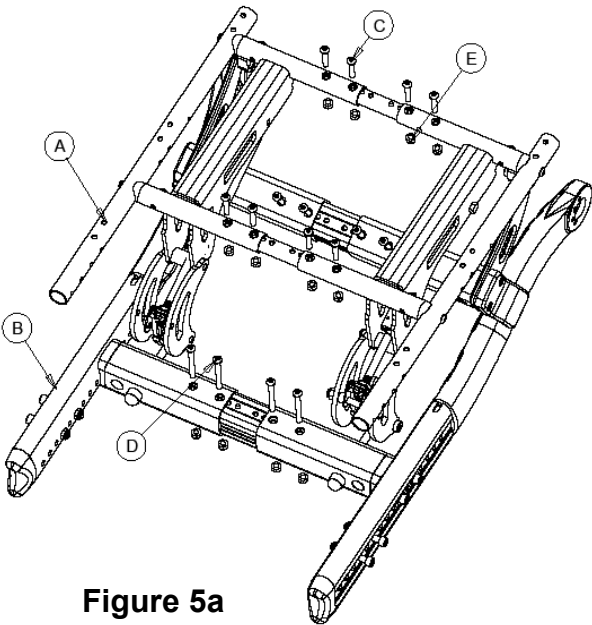


Figure 5a

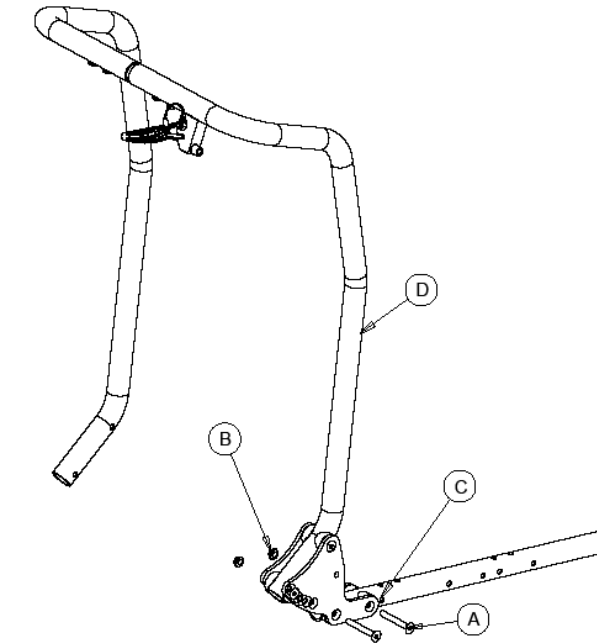


Figure 5d

LOWER CROSSBAR SETUP - CHAIR WIDTH FROM 16" TO 20"

CHAIR WIDTH	16"	17"	18"	19"	20"
SETUP HOLE	1	2	3	4	5


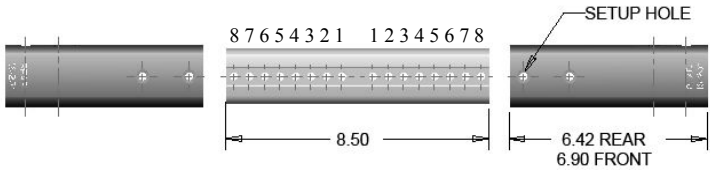


Figure 5b

LOWER CROSSBAR SETUP - CHAIR WIDTH FROM 20" TO 24"

CHAIR WIDTH	20"	21"	22"	23"	24"
SETUP HOLE	1	2	3	4	5

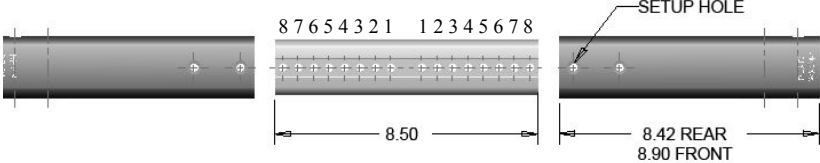

SEATFRAME SETUP - CHAIR WIDTH FROM 16" TO 20"

CHAIR WIDTH	16"	17"	18"	19"	20"
SETUP HOLE	1	2	3	4	5

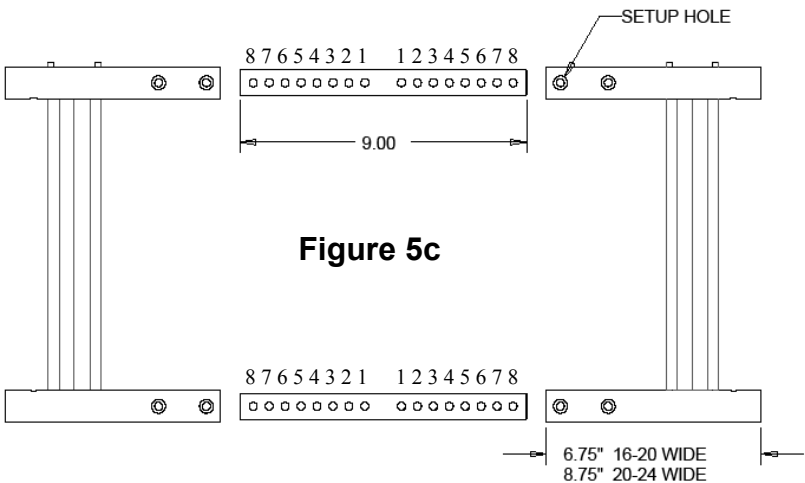


Figure 5c

SEATFRAME SETUP - CHAIR WIDTH FROM 20" TO 24"

CHAIR WIDTH	20"	21"	22"	23"	24"
SETUP HOLE	1	2	3	4	5

Changing Seat Width (Refer to Figure 5a,5b and 5c)

1. Start with either the left side or right side of the wheelchair and remove four (4) 1/4-20 x 1.25" screws('C') and four(4) 1/4-20 nuts('E') from the seat frame('A') (2 in the front and 2 in the back) on the same side.
2. Proceed to the lower frame ('B') and remove four (4) 1/4-20 x 2.25" screws('D') and four (4) 1/4-20 nuts('E') (2 in the front and 2 in the back) on the same side.
3. Grasp the seat frame and lower frame and pull towards oneself observing the holes being uncovered. The Frame can be adjusted by 2" on each side which is 5 holes @ 1/2" center to center on each side. For example to increase the width of the chair by 2" overall each side has to be increased by 1".
4. Place the screws and nuts in their respective holes and hand tighten. The 1/4-20 x 1.25" screws attach to the seat frame while the 1/4-20 x 2.25" screws attach to the lower frame. **See Figure 5b and 5c for hole locations on the crossbar and seatframe.**
5. Now start with the side (either left or right) which remains to be adjusted and remove four (4) 1/4-20 x 1.25" screws('C') and four (4) 1/4-20 nuts('E') from the seat frame('A') as well as remove four (4) 1/4-20 x 2.25" screws('D') and four (4) 1/4-20 nuts ('E') from the lower frame.
6. Repeat step # 3 and #4 and check the overall size of the chair before fully tightening the screws and nuts.

Changing Seat Depth (Refer to Figure 5d)

1. The Capella wheelchair has 3 different depths of seat rail which come with the standard configuration. The depth of the seat rails are 15"-16", 17"-18", and 19"-20". Remove four (4) 1/4-20 x 2" flat head screws('A') and four(4) 1/4-20 nuts('B') from the bottom side of the back cane bracket('C').
2. Move the stroller handle assembly ('D') backwards 1" (for example if one wants to change the depth from 15" to 16") towards the back of the chair and re-insert the four (4) 1/4-20 x 2" flat head screws and four (4) 1/4-20 nuts into the 'new' holes. Tighten the screws and nuts and check the stability of the stroller handle assembly.
3. As the depth of the chair increases or decreases the horizontal location of the wheel must also change. **(Refer Changing Horizontal Location of Rear Wheels).**
4. If one would like to change the depth of the wheelchair from 15"-16" to 19"-20" then the seat rails need to be replaced.

Changing Horizontal Location of Rear Wheels (Refer to Figure 5e)

1. Remove the Rear wheels following the instructions in section "Removing/Installing Rear Wheels".
2. Remove two(2) 5/16-18 x 2" socket cap screws and two(2) 5/16-18 nuts from the wheel lock/anti-tip assembly.
3. If the seat depth is increased by 1" (from 15" to 16") move the wheel lock/anti-tip assembly backwards by 1" and insert the two(2) 5/16-18 x 2" socket cap screws and two(2) 5/16-18 nuts into the 'new' hole locations. If the seat depth is decreased by 1" (from 16" to 15") move the wheel lock/anti -tip assembly forwards by 1" and insert the

screws and nuts.

4. Repeat steps #5 and #6 for the other side of the chair and ensure the screws on both sides of the wheelchair holding the wheel lock/anti-tip assembly are inline.

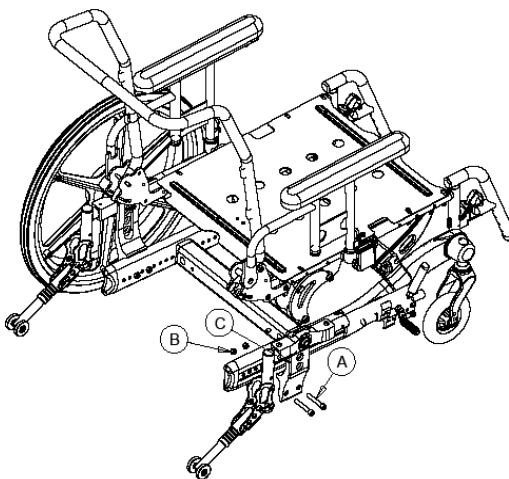


Figure 5e

Changing Vertical Location of Rear Wheels (See section marked SEAT TO FLOOR)

Removing/Installing Rear Wheels(refer to figure 6)

1. Push in the mushroom head ('A') of the quick-release axle ('B') and pull the axle out through the opening in the center of the rear wheel ('C') and shaft spacer ('D').
2. Push in the detent pin of the quick-release axle again and pull the axle out of the rear wheel and shaft spacer.
3. To reinstall the rear wheel onto the wheelchair, reverse steps 1 to 2.
4. If using the Push-to-Lock mechanism ensure the wheel locks properly. See the section ADJUSTING THE WHEEL LOCKS.


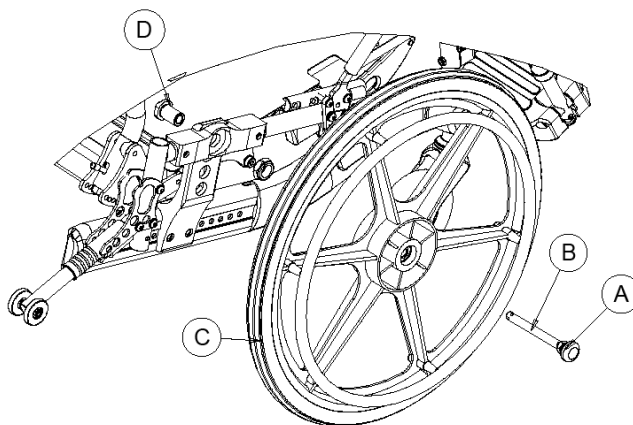


Figure 6

Make sure the detent pin and locking pins of the quick-release axle are fully released BEFORE operating the wheelchair.

The locking pins MUST be protruding past the inside of the rear wheel axle bushing for a positive lock. Keep locking pins clean.

Future Mobility Healthcare Inc. recommends inserting quick-release axles with the detent pin to the inside of the wheelchair to prevent accidental release during contact activities.

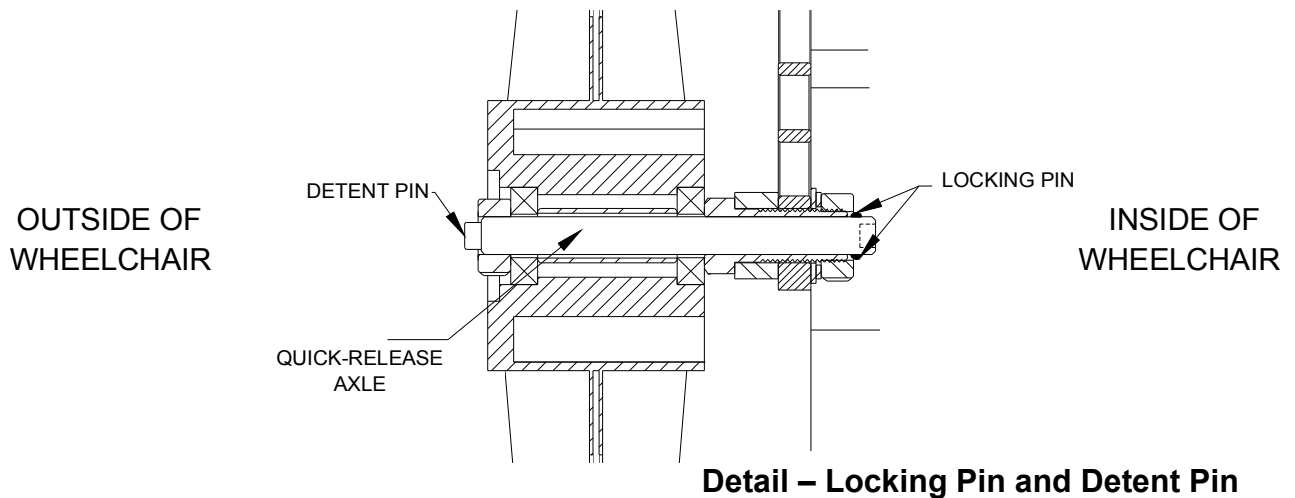


USING THE WHEEL LOCKS(refer to Figure 7a)

WARNING

DO NOT attempt to stop a moving wheelchair with the wheel locks. **WHEEL LOCKS ARE NOT BRAKES** – otherwise injury or damage may occur.

1. Ensure the wheelchair is not moving before engaging the wheel locks.
2. Perform one (1) of the following:
 - a. Push-to-Lock – to engage, push the wheel lock handle ('A') forward
3. Disengage the wheel locks by reversing STEP 2.

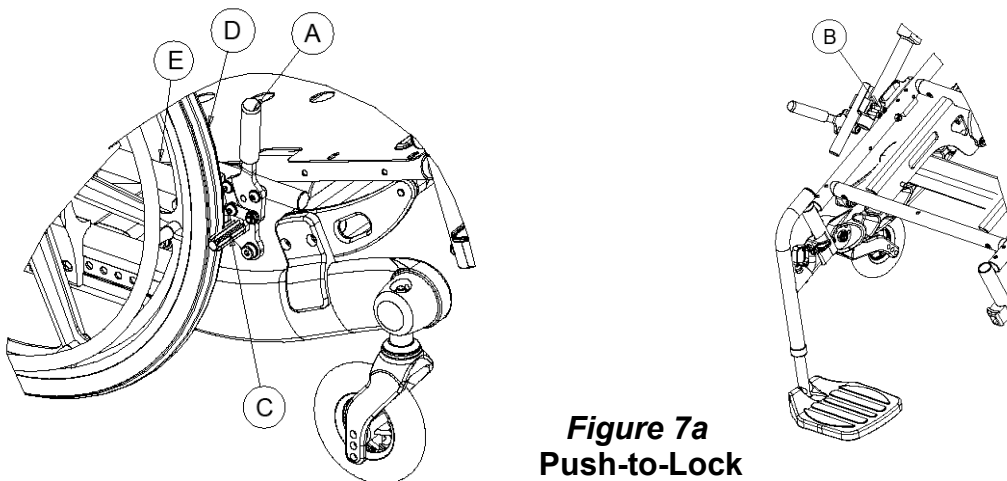

ADJUSTING THE PATIENT OPERATED WHEEL LOCKS

1. Disengage the wheel locks ('A').
2. Loosen the two (2) cap screws ('B') which hold the wheel lock clamps in the above figure that secure the wheel lock to the wheelchair frame.
3. Reposition the wheel lock along the wheel lock tube ('E') so that when engaged, the wheel lock brake shoe ('C') embeds the tire ('D') 1/8" (3/16" for pneumatic tires) and **HOLDS** the occupied wheelchair in place when pushed.
4. Securely tighten the two cap screws securing the wheel lock to the wheelchair frame.
5. Engage the wheel lock.
6. Measure the distance the wheel lock is embedded into the tire
7. Repeat STEPS 1 to 6 until the wheel lock brake shoe embeds the tire and **HOLDS** the occupied wheelchair in place when pushed.
8. Engage both wheel locks and ensure the occupied wheelchair is held in place when pushed.

FRONT CASTERS

WARNING


After **ANY** adjustments, repair or service and **BEFORE** use, make sure all attaching hardware is tightened securely – otherwise injury or damage may occur.

INSTALLING/REPLACING FRONT CASTERS AND FORKS


REPLACING/INSTALLING FRONT CASTERS – refer to Figure 8a

Note: This procedure can be performed if replacing the exact same size front caster.

1. Remove the nut ('A'), bolt ('B') and spacers ('C') which secure the caster to the fork assembly ('D').
2. Remove the front caster.
3. Slide in the new 5,6, or 8 inch front caster into the space between the fork.
4. Tighten the nut, bolt and spacers to secure the caster into place.
5. Repeat this procedure for the other front caster wheel assembly if required.
6. To properly tighten caster journal system and guard against flutter, perform the following check:
 - a. Tip back of wheelchair to floor
 - b. Pivot both forks and casters to top of their arc simultaneously.
 - c. Let casters drop to bottom of arc (wheels should swing once to one-side, then immediately rest in a straight downward position).
 - d. Adjust locknuts according to freedom of caster swing.
7. Test wheelchair for maneuverability.
8. Re-adjust locknuts if necessary and repeat STEPS 6 to 7 until correct.
9. Snap dust cover over the locknut and stem.

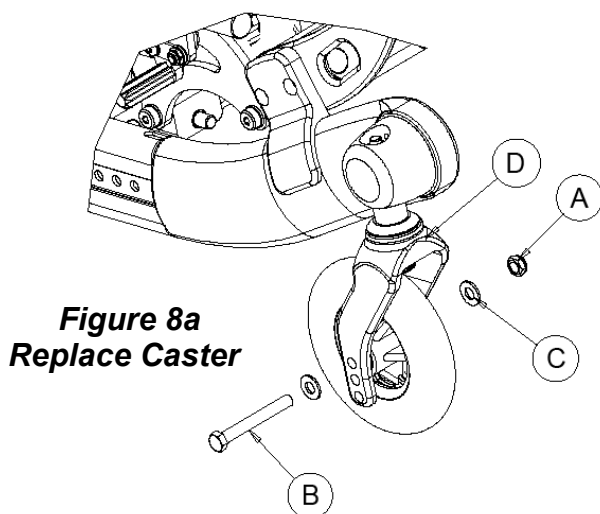

**Figure 8a
Replace Caster**

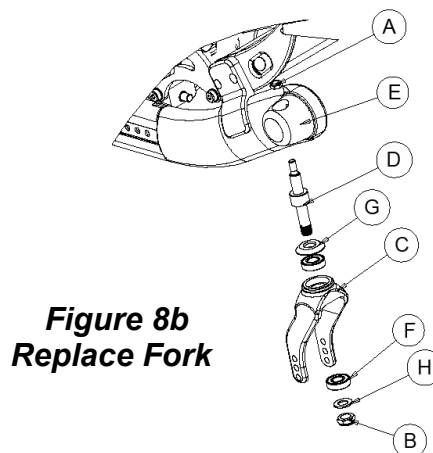

**Figure 8b
Replace Fork**

REPLACING/INSTALLING FRONT FORKS – refer to Figure 8b

1. Remove the locknut ('A') which secures the caster stem bolt ('D') to the caster housing ('E').
2. Drop the front caster and fork assembly out of the caster housing.
3. Remove the locknut ('B') which secures the caster stem to the fork ('C')
4. Also Remove the bearing ('F'), Fork spacer seal ('G') and nylon spacer ('H')
5. Follow REPLACING/INSTALLING FRONT CASTERS guide in this section to remove the front caster from the fork assembly.
6. Place the caster in the new fork and assemble with the caster nut and bolt and spacer.
7. Secure the new fork and caster assembly in place by using the locknut and caster stem bolt.
8. Follow steps 5 to 9 of REPLACING/INSTALLING FRONT CASTERS guide to properly install the new fork assembly and prevent fluttering.

REPLACING/INSTALLING FRONT FORK ASSEMBLY – refer to Figure 8c and 8d

1. Insert the caster fork assembly into the wheelchair main side rail opening as shown in Figure 6c.
2. **Ensure the caster fork is 90° perpendicular to the flat surface**
3. Insert the 1/4-20 x 1" socket head cap screw ('A') through the wheelchair main side rail opening and caster stem housing as shown in Figure 6c. The socket head cap screw is used to hold the fork stem in place so that when the angle of the fork assembly is adjusted it is maintained. The fork angle can be set at 0°, 2°, 4° or 6° as shown in Figure 6d.
4. Also secure the fork assembly to the main side rail by using a 1/2-20 lock nut ('B') as shown in Figure 6c.
5. Install the front caster by following REPLACING/INSTALLING FRONT CASTERS guide In the above section.

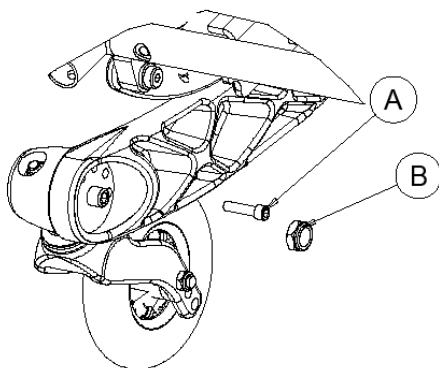


Figure 8c
Replace Fork
Assembly

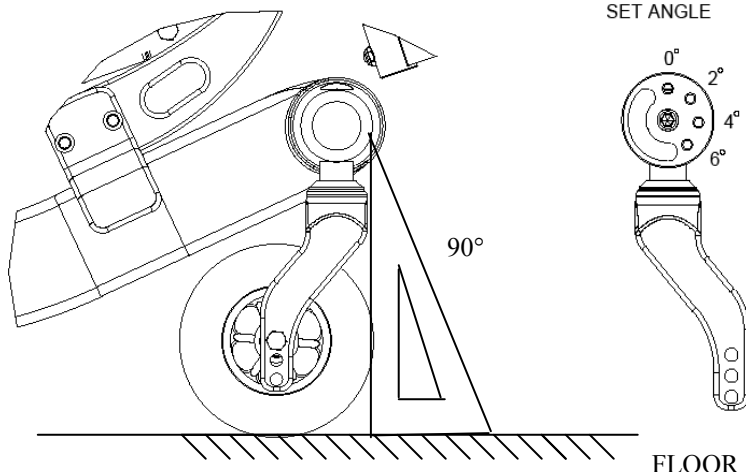



Figure 8d
Set Fork Angle and Ensure 90° to Surface

SEAT TO FLOOR

1. Follow the values in the table below to achieve the desired seat to floor height.
2. Ensure that no weight is on the wheelchair while performing the operations
3. Ensure after changing the seat to floor height that the wheelchair is still operational and that no wheel flutter or instability occurs.
4. Ensure after changing the seat to floor height the caster stem fork is 90° (perpendicular) to the floor.

STF	REAR	AXLE	FORK	CASTER	SPACER	STEM
FRONT	WHEEL	PLATE		SIZE	SIZE	SIZE
(inches)	(inches)	HOLE	HOLE	(inches)	(inches)	
13.00	16	3	1	5	NONE	SHORT
13.00	20	1	1	5	NONE	SHORT
13.50	16	3	2	5	NONE	SHORT
13.50	16	3	1	6	NONE	SHORT
13.50	20	2	2	5	NONE	SHORT
13.50	20	2	1	6	NONE	SHORT
14.00	16	4	3	5	NONE	SHORT
14.00	16	4	2	6	NONE	SHORT
14.00	20	2	3	5	NONE	SHORT
14.00	20	2	2	6	NONE	SHORT
14.00	22	1	3	5	NONE	SHORT
14.00	22	1	2	6	NONE	SHORT
14.50	16	4	3	6	NONE	SHORT
14.50	20	2	3	6	NONE	SHORT
14.50	22	1	3	6	NONE	SHORT
15.00	20	3	3	7	NONE	SHORT
15.00	22	2	3	7	NONE	SHORT
15.00	24	1	3	7	NONE	SHORT

STF	REAR	AXLE	FORK	CASTER	SPACER	STEM
FRONT	WHEEL	PLATE	HOLE	SIZE	SIZE	SIZE
(inches)	(inches)	HOLE	HOLE	(inches)	(inches)	
15.50	20	3	1	5	2.5"	LONG
15.50	22	2	1	5	2.5"	LONG
15.50	24	1	1	5	2.5"	LONG
16.00	20	4	2	5	2.5"	LONG
16.00	20	4	1	6	2.5"	LONG
16.00	22	3	2	5	2.5"	LONG
16.00	22	3	1	6	2.5"	LONG
16.00	24	2	2	5	2.5"	LONG
16.00	24	2	1	6	2.5"	LONG
16.50	20	4	3	5	2.5"	LONG
16.50	20	4	2	6	2.5"	LONG
16.50	22	3	3	5	2.5"	LONG
16.50	22	3	2	6	2.5"	LONG
16.50	24	2	3	5	2.5"	LONG
16.50	24	2	2	6	2.5"	LONG
17.00	22	4	3	6	2.5"	LONG
17.00	22	4	2	7	2.5"	LONG
17.00	24	3	3	6	2.5"	LONG
17.00	24	3	2	7	2.5"	LONG
17.50	22	4	3	7	2.5"	LONG
17.50	24	3	3	7	2.5"	LONG
18.00	24	4	3	8	2.5"	LONG
19.00	24	4	3	8	2.5"	LONG

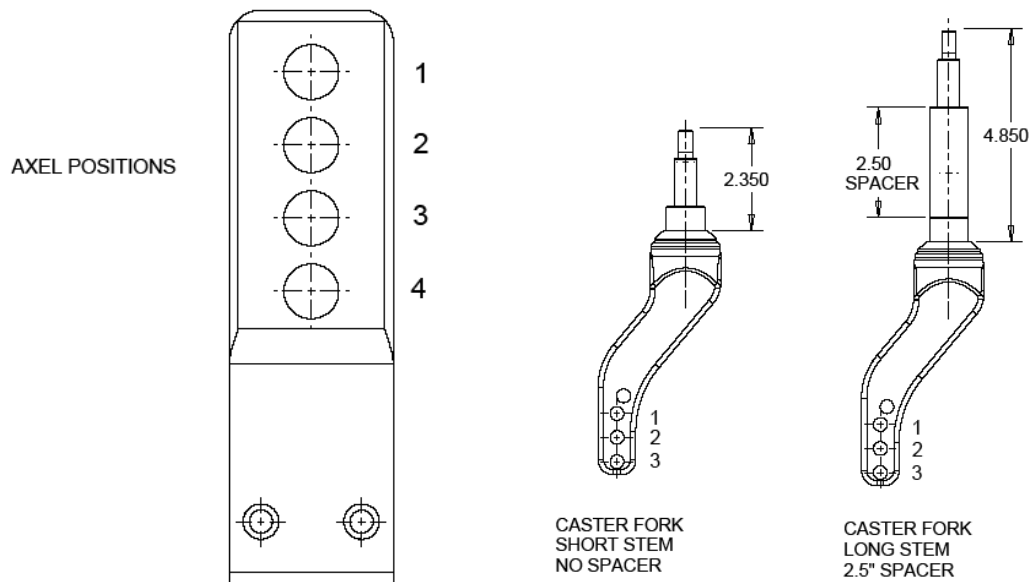

ANTI-TIPPERS/WHEEL LOCKS

WARNING

After **ANY** adjustments, repair or service and **BEFORE** use, make sure all attaching hardware is tightened securely – otherwise injury or damage may occur.

INSTALLING/ADJUSTING THE ANTI-TIPPERS

WARNING

Anti-tippers are specific to the different seat-to floor angles and/or seat-to-floor heights. Refer to the chart in this section of the manual for correct usage and adjustment. If these requirements **CANNOT** be achieved, **DO NOT** use the wheelchair. Contact a qualified technician. If changing the seat-to-floor height with or without a change to seat-to-floor angle, the correct anti-tippers **MUST** be used to maintain a 1 ½" to 2" ground clearance.

Seat-to-floor angle of 3 degrees to 6 degrees: if so equipped, anti-tippers **MUST** be attached at all times. In as much as the anti-tippers are an option on this wheelchair (you may order with or without the anti-tippers), Future Mobility Healthcare Inc. strongly recommends ordering the anti-tippers as a safeguard for the wheelchair user.


Anti-tippers **MUST** be fully engaged and snap buttons fully protruding out of adjustment holes.

Ensure both anti-tippers are adjusted to the same mounting hole.

INSTALLING ANTI-TIPPERS

1. Bolt the anti-tip assembly ('B') onto the anti tip mount tube weldment ('A') at the desired height
2. Measure the distance between the bottom of the anti-tipper wheels and the ground/floor.
3. If the distance between the bottom of the anti-tipper wheels and the ground/floor is not 1-1/2" to 2" adjust the anti tipper location.


*Note: A 1 1/2" to 2" clearance between the bottom of the anti-tipper wheels and the ground/floor **MUST** be maintained at all times.*

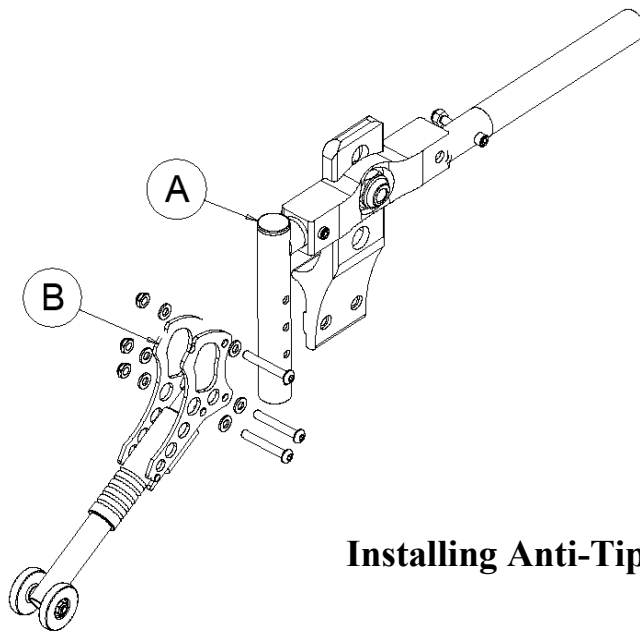

Installing Anti-Tip Assembly

FLIPPING THE ANTI-TIPPERS

1. Pull the sleeve ('A') downwards on the tube to release the pin engagement
2. Rotate the Anti-tip assembly ('B') upwards until the pin engagement is at the top position.

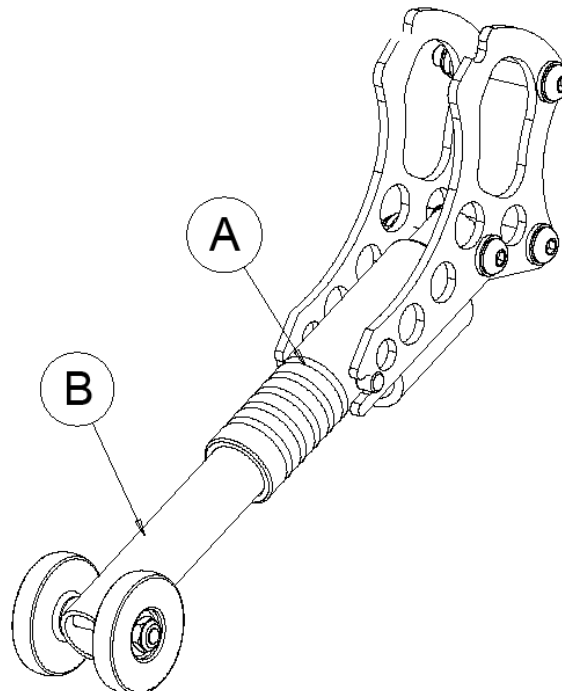

WARRANTY

This warranty is extended only to the original purchaser/user of our products.

Future Mobility Healthcare Inc. warrants its Capella Wheelchair to be free from defects in materials and workmanship for one (1) year from the date of purchase. The frame is warranted for the lifetime of the original purchaser/user. The Orion Back is warranted for three (3) years and the upholstery for 90 days, upon normal usage by original purchaser. If within this warranty period the product shall be proven to be defective, such product shall be repaired or replaced, at Future Mobility Healthcare Inc. discretion. Future Mobility Healthcare Inc. sole obligation and your exclusive remedy under this warranty shall be limited to the repair and/or replacement of the product or its parts. This warranty does not include any labour or shipping charges incurred in replacement part installation or repair of any product.

For warranty service, please contact the dealer from whom you purchased your Future Mobility Healthcare Inc. . In the event you do not receive satisfactory warranty service, please write directly to Future Mobility Healthcare Inc. at 3223 Orlando Drive, Mississauga, Ont. L4V 1C5. Provide the dealer's name, address, model number, date of purchase and indicate the nature of the defect.

DO NOT return products to Future Mobility Healthcare Inc. without our prior consent. The defective unit or parts must be returned for warranty inspection within thirty (30) days of the return authorization date. (Future Mobility Healthcare Inc. will issue a return authorization number). Please prepay all shipping charges; C.O.D. shipments will be refused.

LIMITATIONS and EXCLUSIONS: This warranty shall not apply to problems arising from normal wear or failure to adhere to the enclosed instructions. Products subjected to negligence, accident, improper usage, maintenance or storage; or products modified without Future Mobility Healthcare Inc. written consent including, but not limited to: modification through the use of any unauthorized parts or attachments; products damaged by reason or repairs made to any component without the specific consent of Future Mobility Healthcare Inc., or products repaired by anyone other than a Future Mobility Healthcare Inc. dealer. Such evaluation shall be determined by Future Mobility Healthcare Inc..

The foregoing warranty is exclusive and in lieu of all other expressed warranties. It shall not extend beyond the duration of the expressed warranty provided herein and the remedy for violations of any implied warranty shall be limited to repair or replacement of the defective product pursuant to the terms contained herein. Future Mobility Healthcare Inc. shall not be liable for any consequential or incidental damages whatsoever.

This warranty shall be extended to comply with all state laws and requirements.

Future Mobility Healthcare Inc.
3223 Orlando Drive
Mississauga, ON L4V 1C5
905 671-1661 – 1 888 737-4011
Fax: 905 671-3377
e-mail: customersupport@future-mobility.com
www.future-mobility.com

